

Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.

Regd. No. G-64823


**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44604

**Dr. KETAN DAVDA**

M.B.B.S., M.D.

**PARIKH NURSING HOME**  
(MEDICAL - I.C.C.U.)  
(0260) 2462602

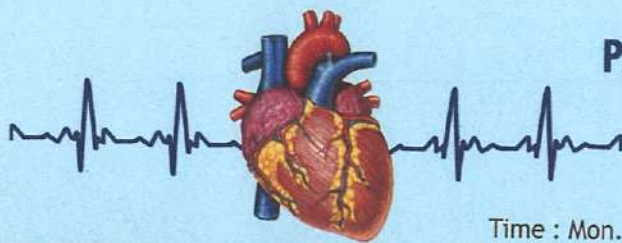

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm  
Sat : 9.00 am - 1.00 pm

<b>Objective</b>	The Objective of this study is to establish the clinical performance by evaluating diagnostic sensitivity and diagnostic specificity of MERISCREEN COVID-19 Antigen Detection kit at Independent site to be used for qualitative detection of SARS-CoV-2 infection in human.		
<b>Scope</b>	The scope of this evaluation is applicable for MERISCREEN COVID-19 Antigen Detection test only.		
<b>Test Name</b>	Meriscreen COVID-19 Antigen Detection Test		
<b>Start Date</b>	9 <sup>th</sup> Dec 2020	<b>End Date</b>	30 <sup>th</sup> Jan 2021

Batch No.	MFG.	EXP.
MI112026	2020/11	2021/10
MI112033	2020/11	2021/10
MI122019	2020/12	2021/11

S.No.	Specimen Type	Specimen ID	PCR Result		COVID-19 Antigen test		
			CT value	Status	C	T	Status
1	Nasopharyngeal	Sample 1	15	Positive	4+	4+	Positive
2	Nasopharyngeal	Sample 2	17	Positive	4+	4+	Positive
3	Nasopharyngeal	Sample 3	26	Positive	4+	2+	Positive
4	Nasopharyngeal	Sample 4	26	Positive	4+	1+	Positive
5	Nasopharyngeal	Sample 5	22	Positive	4+	3+	Positive
6	Nasopharyngeal	Sample 6	26	Positive	4+	2+	Positive
7	Nasopharyngeal	Sample 7	21	Positive	4+	3+	Positive
8	Nasopharyngeal	Sample 8	19	Positive	4+	4+	Positive
9	Nasopharyngeal	Sample 9	28	Positive	4+	2+	Positive
10	Nasopharyngeal	Sample 10	24	Positive	4+	3+	Positive
11	Nasopharyngeal	Sample 11	21	Positive	4+	2+	Positive
12	Nasopharyngeal	Sample 12	26	Positive	4+	1+	Positive
13	Nasopharyngeal	Sample 13	25	Positive	4+	1+	Positive
14	Nasopharyngeal	Sample 14	26	Positive	4+	2+	Positive
15	Nasopharyngeal	Sample 15	25	Positive	4+	2+	Positive
16	Nasopharyngeal	Sample 16	29	Positive	4+	2+	Positive
17	Nasopharyngeal	Sample 17	24	Positive	4+	3+	Positive

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm  
HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064

Regd. No. G-5945


**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.

Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44604

**Dr. KETAN DAVDA**  
M.B.B.S., M.D. (Pathology)

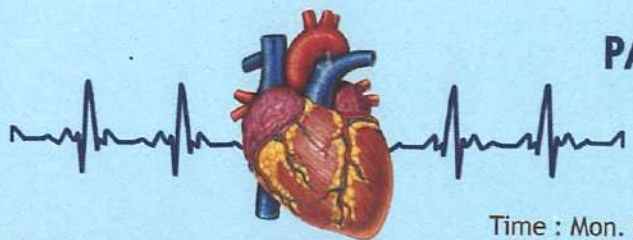

**PARIKH NURSING HOME**  
(MEDICAL - I.C.C.U.)  
(0260) 2462602

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm  
Sat. 9.00 am - 1.00 pm

18	Nasopharyngeal	Sample 18	14	Positive	4+	4+	Positive
19	Nasopharyngeal	Sample 19	26	Positive	4+	1+	Positive
20	Nasopharyngeal	Sample 20	23	Positive	4+	2+	Positive
21	Nasopharyngeal	Sample 21	25	Positive	4+	3+	Positive
22	Nasopharyngeal	Sample 22	27	Positive	4+	1+	Positive
23	Nasopharyngeal	Sample 23	17	Positive	4+	4+	Positive
24	Nasopharyngeal	Sample 24	24	Positive	4+	2+	Positive
25	Nasopharyngeal	Sample 25	15	Positive	4+	4+	Positive
26	Nasopharyngeal	Sample 26	26	Positive	4+	2+	Positive
27	Nasopharyngeal	Sample 27	24	Positive	4+	2+	Positive
28	Nasopharyngeal	Sample 28	24	Positive	4+	3+	Positive
29	Nasopharyngeal	Sample 29	25	Positive	4+	3+	Positive
30	Nasopharyngeal	Sample 30	28	Positive	4+	2+	Positive
31	Nasopharyngeal	Sample 31	16	Positive	4+	4+	Positive
32	Nasopharyngeal	Sample 32	27	Positive	4+	3+	Positive
33	Nasopharyngeal	Sample 33	26	Positive	4+	1+	Positive
34	Nasopharyngeal	Sample 34	22	Positive	4+	2+	Positive
35	Nasopharyngeal	Sample 35	25	Positive	4+	1+	Positive
36	Nasopharyngeal	Sample 36	24	Positive	4+	3+	Positive
37	Nasopharyngeal	Sample 37	26	Positive	4+	1+	Positive
38	Nasopharyngeal	Sample 38	14	Positive	4+	4+	Positive
39	Nasopharyngeal	Sample 39	24	Positive	4+	1+	Positive
40	Nasopharyngeal	Sample 40	32	Positive	4+	+	Positive
41	Nasopharyngeal	Sample 41	17	Positive	4+	4+	Positive
42	Nasopharyngeal	Sample 42	29	Positive	4+	+	Positive
43	Nasopharyngeal	Sample 43	27	Positive	4+	1+	Positive
44	Nasopharyngeal	Sample 44	25	Positive	4+	1+	Positive

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm  
HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064

Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.


Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44604

**Dr. KETAN DAVDA**

M.B.B.S., M.D. Pathology

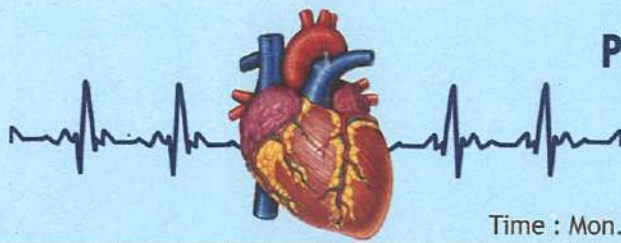

**PARIKH NURSING HOME**  
(MEDICAL - I.C.C.U.)  
(0260) 2462602

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm  
Sat. 8.00 am - 1.00 pm

45	Nasopharyngeal	Sample 45	27	Positive	4+	1+	Positive
46	Nasopharyngeal	Sample 46	29	Positive	4+	+	Positive
47	Nasopharyngeal	Sample 47	26	Positive	4+	2+	Positive
48	Nasopharyngeal	Sample 48	28	Positive	4+	2+	Positive
49	Nasopharyngeal	Sample 49	31	Positive	4+	+	Positive
50	Nasopharyngeal	Sample 50	30	Positive	4+	+	Positive
51	Nasopharyngeal	Sample 51	20	Positive	4+	4+	Positive
52	Nasopharyngeal	Sample 52	27	Positive	4+	+	Positive
53	Nasopharyngeal	Sample 53	29	Positive	4+	+	Positive
54	Nasopharyngeal	Sample 54	25	Positive	4+	+	Positive
55	Nasopharyngeal	Sample 55	22	Positive	4+	3+	Positive
56	Nasopharyngeal	Sample 56	26	Positive	4+	1+	Positive
57	Nasopharyngeal	Sample 57	25	Positive	4+	1+	Positive
58	Nasopharyngeal	Sample 58	30	Positive	4+	+	Positive
59	Nasopharyngeal	Sample 59	29	Positive	4+	+	Positive
60	Nasopharyngeal	Sample 60	24	Positive	4+	2+	Positive
61	Nasopharyngeal	Sample 61	19	Positive	4+	4+	Positive
62	Nasopharyngeal	Sample 62	30	Positive	4+	+	Positive
63	Nasopharyngeal	Sample 63	31	Positive	4+	+	Positive
64	Nasopharyngeal	Sample 64	29	Positive	4+	+	Positive
65	Nasopharyngeal	Sample 65	27	Positive	4+	+	Positive
66	Nasopharyngeal	Sample 66	22	Positive	4+	4+	Positive
67	Nasopharyngeal	Sample 67	26	Positive	4+	1+	Positive
68	Nasopharyngeal	Sample 68	29	Positive	4+	+	Positive
69	Nasopharyngeal	Sample 69	27	Positive	4+	+	Positive
70	Nasopharyngeal	Sample 70	26	Positive	4+	+	Positive
71	Nasopharyngeal	Sample 71	29	Positive	4+	+	Positive

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm

HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064

Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.


Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44604

**Dr. KETAN DAVDA**

M.B.B.S., M.D. (Pathology)

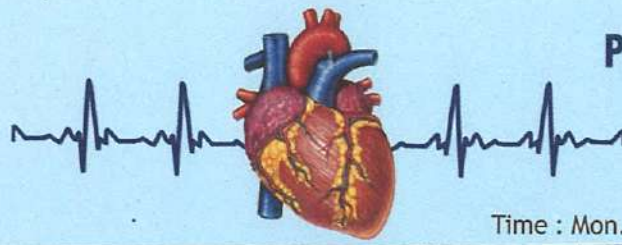

**PARIKH NURSING HOME**  
(MEDICAL - I.C.C.U.)  
(0260) 2462602

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm  
9.00 am - 1.00 pm

72	Nasopharyngeal	Sample 72	29	Positive	4+	+	Positive
73	Nasopharyngeal	Sample 73	28	Positive	4+	+	Positive
74	Nasopharyngeal	Sample 74	27	Positive	4+	+	Positive
75	Nasopharyngeal	Sample 75	25	Positive	4+	2+	Positive
76	Nasopharyngeal	Sample 76	26	Positive	4+	1+	Positive
77	Nasopharyngeal	Sample 77	23	Positive	4+	2+	Positive
78	Nasopharyngeal	Sample 78	28	Positive	4+	+	Positive
79	Nasopharyngeal	Sample 79	26	Positive	4+	1+	Positive
80	Nasopharyngeal	Sample 80	27	Positive	4+	1+	Positive
81	Nasopharyngeal	Sample 81	25	Positive	4+	2+	Positive
82	Nasopharyngeal	Sample 82	34	Positive	4+	0	Negative
83	Nasopharyngeal	Sample 83	28	Positive	4+	+	Positive
84	Nasopharyngeal	Sample 84	26	Positive	4+	+	Positive
85	Nasopharyngeal	Sample 85	27	Positive	4+	+	Positive
86	Nasopharyngeal	Sample 86	26	Positive	4+	+	Positive
87	Nasopharyngeal	Sample 87	33	Positive	4+	0	Negative
88	Nasopharyngeal	Sample 88	24	Positive	4+	2+	Positive
89	Nasopharyngeal	Sample 89	27	Positive	4+	1+	Positive
90	Nasopharyngeal	Sample 90	27	Positive	4+	+	Positive
91	Nasopharyngeal	Sample 91	23	Positive	4+	1+	Positive
92	Nasopharyngeal	Sample 92	33	Positive	4+	0	Negative
93	Nasopharyngeal	Sample 93	25	Positive	4+	1+	Positive
94	Nasopharyngeal	Sample 94	24	Positive	4+	2+	Positive
95	Nasopharyngeal	Sample 95	26	Positive	4+	2+	Positive
96	Nasopharyngeal	Sample 96	27	Positive	4+	+	Positive
97	Nasopharyngeal	Sample 97	23	Positive	4+	1+	Positive
98	Nasopharyngeal	Sample 98	25	Positive	4+	1+	Positive

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm

HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064

Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.


Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44604

**Dr. KETAN DAVDA**

M.B.B.S., M.D. (Pathology)

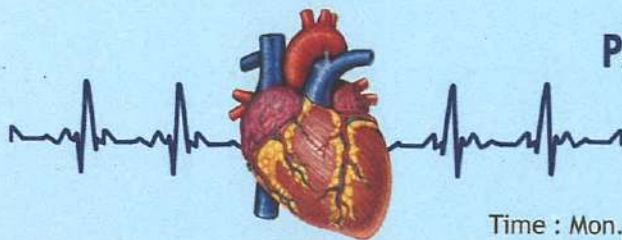

**PARIKH NURSING HOME**  
(MEDICAL - I.C.C.U.)  
(0260) 2462602

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm  
Sat. 9.00 am - 1.00 pm

99	Nasopharyngeal	Sample 99	26	Positive	4+	1+	Positive
100	Nasopharyngeal	Sample 100	29	Positive	4+	+	Positive
101	Nasopharyngeal	Sample 101	28	Positive	4+	1+	Positive
102	Nasopharyngeal	Sample 102	27	Positive	4+	+	Positive
103	Nasopharyngeal	Sample 103	26	Positive	4+	+	Positive
104	Nasopharyngeal	Sample 104	29	Positive	4+	+	Positive
<b>Total No. of PCR Positive Samples</b>		<b>104</b>					
<b>No. of sample Positive in Antigen test</b>		<b>101</b>					
<b>Negative sample</b>		<b>03</b>					
<b>Diagnostic Sensitivity</b>		<b>97.20%</b>					

**Diagnostic Specificity:**

S.No.	Specimen Type	Specimen ID	PCR Result	COVID-19 Antigen test		
			Status	C	T	Status
1	Nasopharyngeal	Sample 1	Negative	4+	0	Negative
2	Nasopharyngeal	Sample 2	Negative	4+	0	Negative
3	Nasopharyngeal	Sample 3	Negative	4+	0	Negative
4	Nasopharyngeal	Sample 4	Negative	4+	0	Negative
5	Nasopharyngeal	Sample 5	Negative	4+	0	Negative
6	Nasopharyngeal	Sample 6	Negative	4+	0	Negative
7	Nasopharyngeal	Sample 7	Negative	4+	0	Negative
8	Nasopharyngeal	Sample 8	Negative	4+	0	Negative
9	Nasopharyngeal	Sample 9	Negative	4+	0	Negative
10	Nasopharyngeal	Sample 10	Negative	4+	0	Negative
11	Nasopharyngeal	Sample 11	Negative	4+	0	Negative
12	Nasopharyngeal	Sample 12	Negative	4+	0	Negative

**Nucleus Clinical Laboratory**

**Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm**  
**HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064**

Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.


Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44604

**Dr. KETAN DAVDA**

M.B.B.S., M.D. (Pathology)

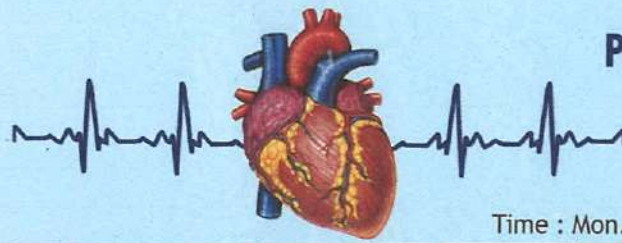

**PARIKH NURSING HOME**  
(MEDICAL - I.C.C.U.)  
(0260) 2462602

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm  
Sat. 9.00 am - 1.00 pm

13	Nasopharyngeal	Sample 13	Negative	4+	0	Negative
14	Nasopharyngeal	Sample 14	Negative	4+	0	Negative
15	Nasopharyngeal	Sample 15	Negative	4+	0	Negative
16	Nasopharyngeal	Sample 16	Negative	4+	0	Negative
17	Nasopharyngeal	Sample 17	Negative	4+	0	Negative
18	Nasopharyngeal	Sample 18	Negative	4+	0	Negative
19	Nasopharyngeal	Sample 19	Negative	4+	0	Negative
20	Nasopharyngeal	Sample 20	Negative	4+	0	Negative
21	Nasopharyngeal	Sample 21	Negative	4+	0	Negative
22	Nasopharyngeal	Sample 22	Negative	4+	0	Negative
23	Nasopharyngeal	Sample 23	Negative	4+	0	Negative
24	Nasopharyngeal	Sample 24	Negative	4+	0	Negative
25	Nasopharyngeal	Sample 25	Negative	4+	0	Negative
26	Nasopharyngeal	Sample 26	Negative	4+	0	Negative
27	Nasopharyngeal	Sample 27	Negative	4+	0	Negative
28	Nasopharyngeal	Sample 28	Negative	4+	0	Negative
29	Nasopharyngeal	Sample 29	Negative	4+	0	Negative
30	Nasopharyngeal	Sample 30	Negative	4+	0	Negative
31	Nasopharyngeal	Sample 31	Negative	4+	0	Negative
32	Nasopharyngeal	Sample 32	Negative	4+	0	Negative
33	Nasopharyngeal	Sample 33	Negative	4+	0	Negative
34	Nasopharyngeal	Sample 34	Negative	4+	0	Negative
35	Nasopharyngeal	Sample 35	Negative	4+	0	Negative
36	Nasopharyngeal	Sample 36	Negative	4+	0	Negative
37	Nasopharyngeal	Sample 37	Negative	4+	0	Negative
38	Nasopharyngeal	Sample 38	Negative	4+	0	Negative
39	Nasopharyngeal	Sample 39	Negative	4+	0	Negative

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm

HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064

Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.


Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44804

**Dr. KETAN DAVDA**

M.B.B.S., M.D. (Pathology)

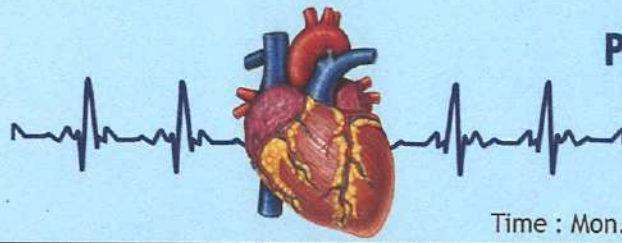

**PARIKH NURSING HOME**  
(MEDICAL - I.C.C.U.)  
(0260) 2462602

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm  
Sat. 8.00 am - 1.00 pm

40	Nasopharyngeal	Sample 40	Negative	4+	0	Negative
41	Nasopharyngeal	Sample 41	Negative	4+	0	Negative
42	Nasopharyngeal	Sample 42	Negative	4+	0	Negative
43	Nasopharyngeal	Sample 43	Negative	4+	0	Negative
44	Nasopharyngeal	Sample 44	Negative	4+	0	Negative
45	Nasopharyngeal	Sample 45	Negative	4+	0	Negative
46	Nasopharyngeal	Sample 46	Negative	4+	0	Negative
47	Nasopharyngeal	Sample 47	Negative	4+	0	Negative
48	Nasopharyngeal	Sample 48	Negative	4+	0	Negative
49	Nasopharyngeal	Sample 49	Negative	4+	0	Negative
50	Nasopharyngeal	Sample 50	Negative	4+	0	Negative
51	Nasopharyngeal	Sample 51	Negative	4+	0	Negative
52	Nasopharyngeal	Sample 52	Negative	4+	0	Negative
53	Nasopharyngeal	Sample 53	Negative	4+	0	Negative
54	Nasopharyngeal	Sample 54	Negative	4+	0	Negative
55	Nasopharyngeal	Sample 55	Negative	4+	0	Negative
56	Nasopharyngeal	Sample 56	Negative	4+	0	Negative
57	Nasopharyngeal	Sample 57	Negative	4+	0	Negative
58	Nasopharyngeal	Sample 58	Negative	4+	0	Negative
59	Nasopharyngeal	Sample 59	Negative	4+	0	Negative
60	Nasopharyngeal	Sample 60	Negative	4+	0	Negative
61	Nasopharyngeal	Sample 61	Negative	4+	0	Negative
62	Nasopharyngeal	Sample 62	Negative	4+	0	Negative
63	Nasopharyngeal	Sample 63	Negative	4+	0	Negative
64	Nasopharyngeal	Sample 64	Negative	4+	0	Negative
65	Nasopharyngeal	Sample 65	Negative	4+	0	Negative
66	Nasopharyngeal	Sample 66	Negative	4+	0	Negative

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm  
HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064

Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.


Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)

Regd. No. G-44604

**Dr. KETAN DAVDA**

M.B.B.S., M.D. (Pathology)

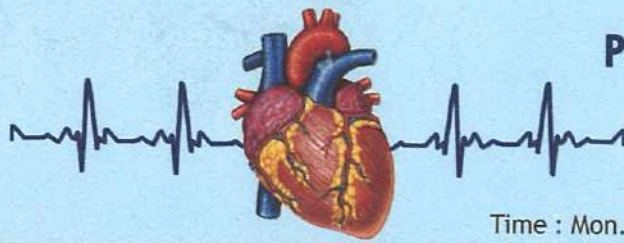

**PARIKH NURSING HOME**

(MEDICAL - I.C.C.U.)

(0260) 2462602

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm

Sat. 9.00 am - 1.00 pm

67	Nasopharyngeal	Sample 67	Negative	4+	0	Negative
68	Nasopharyngeal	Sample 68	Negative	4+	0	Negative
69	Nasopharyngeal	Sample 69	Negative	4+	0	Negative
70	Nasopharyngeal	Sample 70	Negative	4+	0	Negative
71	Nasopharyngeal	Sample 71	Negative	4+	0	Negative
72	Nasopharyngeal	Sample 72	Negative	4+	0	Negative
73	Nasopharyngeal	Sample 73	Negative	4+	0	Negative
74	Nasopharyngeal	Sample 74	Negative	4+	0	Negative
75	Nasopharyngeal	Sample 75	Negative	4+	0	Negative
76	Nasopharyngeal	Sample 76	Negative	4+	0	Negative
77	Nasopharyngeal	Sample 77	Negative	4+	0	Negative
78	Nasopharyngeal	Sample 78	Negative	4+	0	Negative
79	Nasopharyngeal	Sample 79	Negative	4+	0	Negative
80	Nasopharyngeal	Sample 80	Negative	4+	0	Negative
81	Nasopharyngeal	Sample 81	Negative	4+	0	Negative
82	Nasopharyngeal	Sample 82	Negative	4+	0	Negative
83	Nasopharyngeal	Sample 83	Negative	4+	0	Negative
84	Nasopharyngeal	Sample 84	Negative	4+	0	Negative
85	Nasopharyngeal	Sample 85	Negative	4+	0	Negative
86	Nasopharyngeal	Sample 86	Negative	4+	0	Negative
87	Nasopharyngeal	Sample 87	Negative	4+	0	Negative
88	Nasopharyngeal	Sample 88	Negative	4+	0	Negative
89	Nasopharyngeal	Sample 89	Negative	4+	0	Negative
90	Nasopharyngeal	Sample 90	Negative	4+	0	Negative
91	Nasopharyngeal	Sample 91	Negative	4+	0	Negative
92	Nasopharyngeal	Sample 92	Negative	4+	0	Negative
93	Nasopharyngeal	Sample 93	Negative	4+	0	Negative

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm

HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064



Regd. No. G-5945

**Dr. KIRIT N. PARIKH**  
CONSULTANT PHYSICIAN M.D., A.F.I.H.

Regd. No. G-64823

**Dr. RUTA K. PARIKH**  
Consultant in Critical Care M.D. (Anae.)


Regd. No. G-44604

**Dr. KETAN DAVDA**

**PARIKH NURSING HOME**

(MEDICAL - I.C.C.U.)

(0260) 2462602

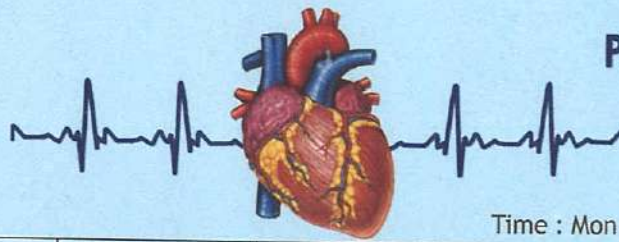

Time : Mon. to Fri. 9.00 am. & 1.00 pm - 7.30 pm


Sat. 9.00 am - 1.00 pm

Sample No.	Sample Type	Sample ID	Result	PCR	Antigen	Final Result
94	Nasopharyngeal	Sample 94	Negative	4+	0	Negative
95	Nasopharyngeal	Sample 95	Negative	4+	0	Negative
96	Nasopharyngeal	Sample 96	Negative	4+	0	Negative
97	Nasopharyngeal	Sample 97	Negative	4+	0	Negative
98	Nasopharyngeal	Sample 98	Negative	4+	0	Negative
99	Nasopharyngeal	Sample 99	Negative	4+	0	Negative
100	Nasopharyngeal	Sample 100	Negative	4+	0	Negative

**Conclusion:**

Diagnostic Sensitivity		Diagnostic Specificity	
PCR Positive	104	PCR Negative	100
Antigen Positive	101	Antigen Positive	0
Antigen Negative	03	Antigen Negative	100
Total No. Of samples	104	Total No. of samples	100
Sensitivity (%)	97.20%	Specificity (%)	100%

Name and Address of Testing Site	Nucleus clinical laboratory PARIKH NURSING HOME, VAPI. ☎ : 0260-2462602
Evaluator's Name and Signature with seal/stamp	DR. K. N. PARIKH <i>[Signature]</i> Dr. Kirit N. Parikh Regd. No. G-5945 M.D.A.F.I.H Parikh Nursing Home Daman Road, Chala-Vapi
Date:	16-2-21

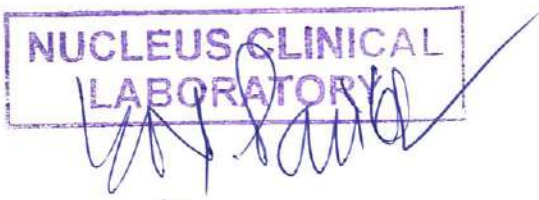

**Nucleus Clinical Laboratory**

Timing : Mon. to Sat. - 8:00 am to 10:00 pm / Sun - 8.00 am to 9.00 pm  
HOME VISIT ON REQUEST FOR NUCLEUS CLINICAL LAB : +91 96876 54064