




Diagnostics

MeriScreen™ | HIV 1-2 WB


One Step qualitative immunochromatographic assay for detection of antibodies against HIV - 1 and HIV - 2 in human serum, plasma and whole blood (Fingerstick & Venous)

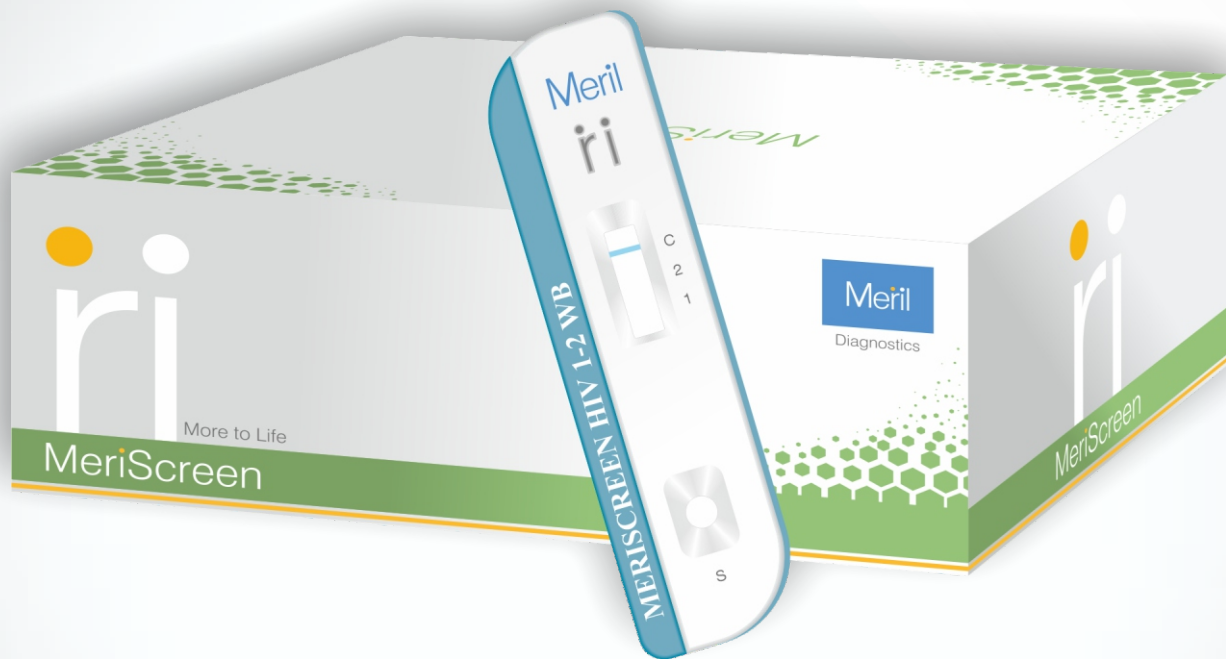

SPECIFICATION

- Immunochromatography assay
- Detection of IgM, IgG and IgA Antibodies
- Specimen Type: Serum, Plasma, Whole blood (Fingerstick & Venous)
- Specimen Volume: 10 μ l Serum/ Plasma or 20 μ l Whole Blood
- Buffer Volume: 3 drops
- Results after 20 minutes but no later than 30 minute
- Diagnostic Sensitivity: 100 %
- Diagnostic Specificity: 100 %
- Shelf life of 24 months

TEST PROTOCOL

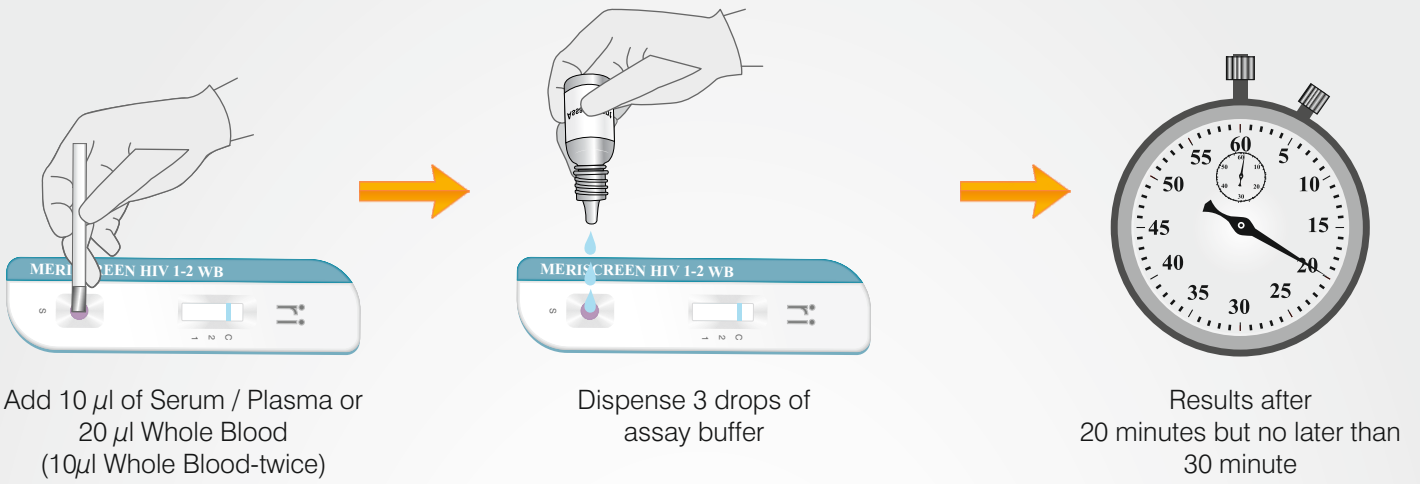

INTERPRETATION

Non-Reactive 	Reactive 	Invalid 			
			Mixed Infection	HIV-1 Infection	HIV-2 Infection

No Pink- Purple Line At 'C' Region Even After Addition of Sample & Assay Buffer Indicates Invalid Test

ORDERING INFORMATION

Product Code	Material Description	Pack Size
HVWRPD-01	MERISCREEN HIV 1-2 WB	30 Test
HVWRPD-02	MERISCREEN HIV 1-2 WB	40 Test



Meril Diagnostics Pvt. Ltd.

D1- D3, Second Floor, Meril Park, Survey No. 135/2/B & 174/2, Muktanand Marg, Chala, Vapi - 396191, India

Tel: +91 260 2408000 | Fax: +91 260 2408025 | Email: askinfo@merillife | mahendra.s@merillife.com