

MeriScreen™

MERISCREEN Malaria Pf HRP-II Ag

Evaluated
& Pre Qualified
by WHO

Meril
Diagnostics



SPECIFICATION


Type	Rapid immunochromatographic test kit
Detected Species	<i>P. falciparum</i>
Specimen	Whole Blood
Specimen transfer device	Inverted Cup
Read Time	20 minutes
Storage Temperature	1-40°C
Shelf Life	24 months
Limit of Detection (LOD)	≥50 parasites/μl

MERISCREEN Malaria Pf HRP-II Ag




Diagnostics

TEST PROTOCOL

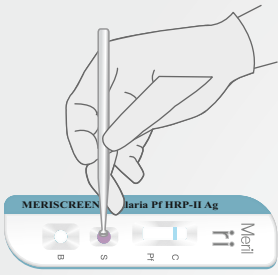

Add 5 μ l of whole blood using Specimen transfer device

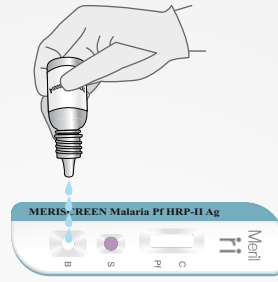

Dispense 4 Drops of Assay buffer



Results in 20 minutes

INTERPRETATION

Positive for <i>P. falciparum</i>	Positive for <i>P. falciparum</i>	Negative	Invalid. Take a new device and repeat the test	Invalid. Take a new device and repeat the test	Invalid. Take a new device and repeat the test	Invalid. Take a new device and repeat the test

ORDERING INFORMATION

PRODUCT CODE	DESCRIPTION	Pack Size
RPWMPH-02	MERISCREEN Malaria Pf HRP-II Ag	10 Tests
RPWMPH-03	MERISCREEN Malaria Pf HRP-II Ag	25 Tests
RPWMPH-04	MERISCREEN Malaria Pf HRP-II Ag	50 Tests

An ISO 13485 & GMP Certified Company



Meril Diagnostics Pvt. Ltd.

D1- D3, Second Floor, Meril Park, Survey No. 135/2/B & 174/2, Muktanand Marg, Chala, Vapi - 396191, India

+91 9624025718 | +91 9429632124 | Email: mahendra.s@merillife.com | askinfo@merillife.com